

MODULAR TABLE SYSTEM



Spinal Surgery Top



Radiolucent Imaging Top



Orthopedic Trauma Top



Advanced Control Base

MIZUHOSI[®]



BRINGING INNOVATION TO THE TABLE.[®]

ortho trauma

FOR COMPLETE ORTHOPEDIC TRAUMA

With the Mizuho OSI Advanced Control Modular Table Base, the Orthopedic Trauma Top provides a complete, powered, orthopedic table system for easy patient positioning and provides for skin or skeletal traction during key trauma procedures.

- HIP PINNINGS
- FEMUR, SUPINE OR LATERAL IM NAILING
- ENDER PINNING
- TIBIA IM NAILING
- HUMERUS IM NAILING

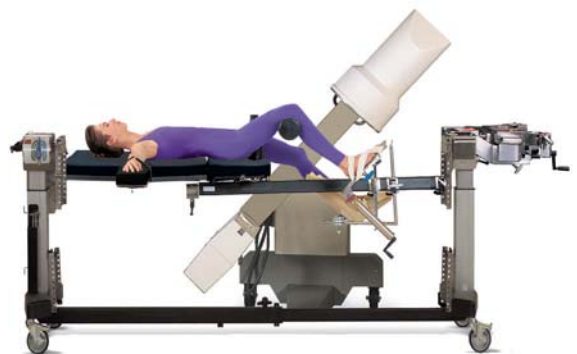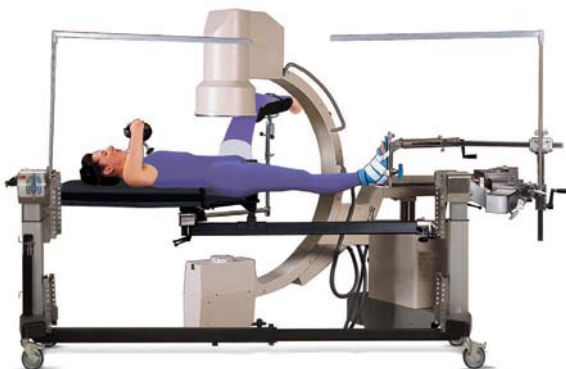

PROCEDURES.



Left Hip Procedure, Bilateral Skin Traction



Right Femoral Nailing, Unilateral Skeletal Traction



Left Femoral Nailing, Skeletal Traction



Left Femoral Nailing, Skin Traction

MAXIMUM RADIOLUCENCY

- Unique composite top enhances radiography
- Allows upper and lower extremity imaging at any angle

ANATOMICALLY CORRECT TRACTION

- Traction Arc provides constant traction on patient's leg for adduction/abduction
- Eliminates need for leg spar on operative side
- Permits unrestricted access for surgeon and C-Arm

DUAL COLUMNS ENHANCE STABILITY

- Patient can be centered for added stability and C-Arm accessibility

LIGHTWEIGHT CONSTRUCTION

- Carbon composite construction adds strength and reduces weight
- Steer casters aid table transport

POWER CONTROLS

- Hand controller facilitates height, trandelenburg / reverse trandelenburg adjustments

ACCESSORY CART

- Specially designed for orthopedic accessories storage and retrieval
- Easy transfer of Traction Arc onto the table



MODULAR BASE TECHNOLOGY

- Accepts Orthopedic Trauma, Spinal, and Imaging Tops
- Creates totally comprehensive orthopedic trauma table capabilities
- Allows supine and prone procedures without removing patient

ORDERING INFORMATION

REF

- 5855 — Orthopedic Trauma Top w/ 500 lbs. Patient Load
- 5803 — Advanced Control Base 120V, 60Hz
- 5803i — Advanced Control Base 220/230V, 50/60Hz
- 5803j — Advanced Control Base 100V, 50/60Hz

SPECIFICATIONS

ORTHOPEDIC TRAUMA TOP

- Table Top Length84 in. (213 cm)
- Table Top Width.....21-1/2 in. (55 cm)
- Patient Weight Capacity500 lbs. (227 kg)

ADVANCED CONTROL BASE

- Table Top to Floor Height.....34 in. – 50 in. (86 cm – 127 cm)
- Base Width32 in. (81 cm)
- Base Length..... 102 in. (259 cm)
- Retractable Base Length64 in. – 102 in. (163 cm – 259 cm)
- Trendelenburg / Reverse Trendelenburg..... 10°
- Lateral Tilt.....25°
- Patient Weight Capacity500 lbs. (227 kg)

OPTIONAL ACCESSORIES

REF

- 5356.....Pivoting Arm Board w/ Tempur-Pedic® Medical Pad
- 5857.....Cross Arm Support
- 5393.....Clark Socket
- 6977-959.....Universal Side Rail Adaptor, 6 in. (15 cm)
- 5858.....Top Patient Isolation Drapes, 12/case
- 5855-865.....Traction Boots, extra-large, 1 pair
- 5855-17.....Traction Boot Liners, large, 24 pack
- 5855-18.....Traction Boot Liners, small, 24 pack
- 5364.....Lateral Arm Board Set
- 5859.....Orthopedic Equipment Cart
- 10203.....OrangeAid® Gel Closed Head Ring, adult
- 10302.....OrangeAid® Gel Armboard Pad, long, extra-wide
- 10503.....OrangeAid® Gel Knee Crutch Pads, medium, set of 2
- 10704.....OrangeAid® Gel Universal Perineal Post Pad
- 10706.....OrangeAid® Gel Universal Traction Boot Liner
- 10920.....OrangeAid® Gel Patient Safety Strap

STANDARD EQUIPMENT

- T Pin
- 90 Degree Pin & Wire Holder
- Traction Slide Assembly
- Tibia Counter Traction Assembly
- Latreal Perineal Post Assembly
- Traction Boot, large
- Traction Boot, small
- Well Leg Board Assembly
- Traction Extension, long, 38 in. (96.52 cm)
- Ilium Post Assembly
- Pad - Counter Traction
- Perineal Post w/Pad & Lock Pin, adult
- Perineal Post w/Pad, child
- Cover, Orthopedic Trauma Top w/o Arc
- Center Support Board
- Traction Unit
- Tibia Traction Upright
- Foot Plate Assembly
- Pad - Lateral Extension
- G.M. Board w/Locking Knob, right
- G.M. Board w/Locking Knob, left
- Safety Strap
- Traction Extension, short
- Accessory Clamp
- Drapery Rod Assembly
- Head End Drapery Post, Weldment
- Foot End Drapery Post, Weldment
- Pad - Lateral Post, lower
- Leg Spar Upright
- Traction Unit Upright, left
- Traction Unit Upright, right
- Radiolucent Leg Spar
- Well Leg Lower Swing Arm
- Lower Leg Support Assembly
- Chart Set Assembly, MTS 5855
- Universal Side Rail Adaptor, 6 in. (15 cm)
- Well Leg Support Assembly
- Traction Boot w/Foot Plate, small, left & right
- Traction Boot w/Foot Plate, large, left & right
- Table Top Assembly
- Transfer Board Assembly
- Tibia Well Leg Support Assembly
- Patient Isolation Drapes, box of 3
- Kirchner Bow/Holder
- Orthopedic Traction Arc Cart
- Tabletop Support Assembly


MIZUHOSI®

 MIZUHO OSI
30031 AHERN AVENUE
UNION CITY, CA 94587-1234/USA
TELEPHONE: 510 - 429 - 1500
INTERNATIONAL: + 1 - 510 - 429 - 1500
TOLL FREE: 800 - 777 - 4674
FAX: 510 - 429 - 8500
WWW.MIZUHOSI.COM
WWW.NEWHIPNEWS.COM



EC REP

EMERGO EUROPE
MOLENSTRAAT 15
2513 BH THE HAGUE
THE NETHERLANDS
TEL: (31) (0) 70 345-5870
FAX: (31) (0) 70 346-7299

Mizuho OSI's ISO 13485 Quality Management System is certified by TÜV Rheinland® LGA Products GmbH.  TÜVRheinland®
Note: Mizuho OSI is constantly improving its products, all specifications are subject to change without notice.
Tempur-Pedic® is a registered trademark of Tempur-Pedic North America, LLC.
Mizuho OSI is a Delaware Corporation. Manufactured in the USA.